


Model 03-206-05-CM

Pro-MAM BR Mammography Phototimer Test Plates - 8 cm

Part Number: 99-8992

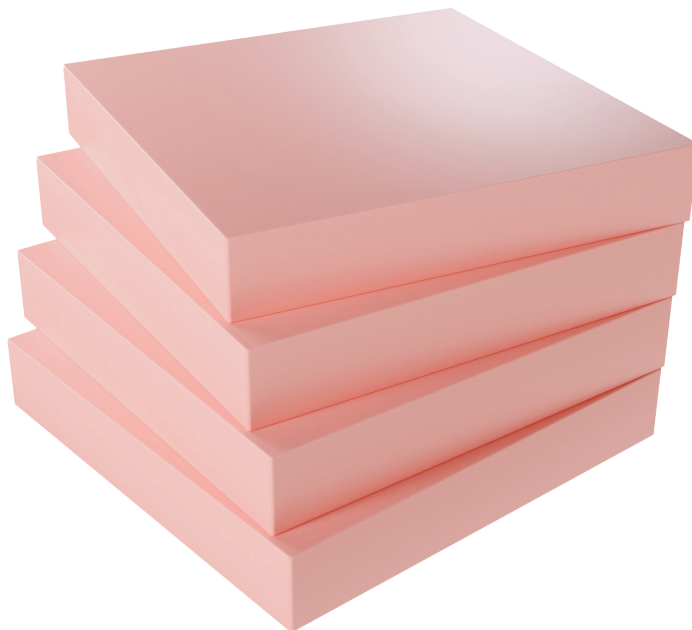

Introduction

The Model 03-206-05-CM Pro-MAM BR Mammography Phototimer Test Plates are designed for accurate assessment of AEC system performance in accordance with American College of Radiology and MQSA recommendations. Unlike acrylic, these testing slabs are manufactured with very tight thickness tolerances and more accurately simulate real breast tissue over the range of energies used in mammography. The Pro-MAM BR (50% glandular / 50% adipose) is most commonly used.

Specifications

- Set of 4 adipose / gland tissues slabs
 - 4 pcs. - 12 x 12 x 2 cm
- Complies with IEC 61223-3-2
- CE certified